
# PROFI CARE®

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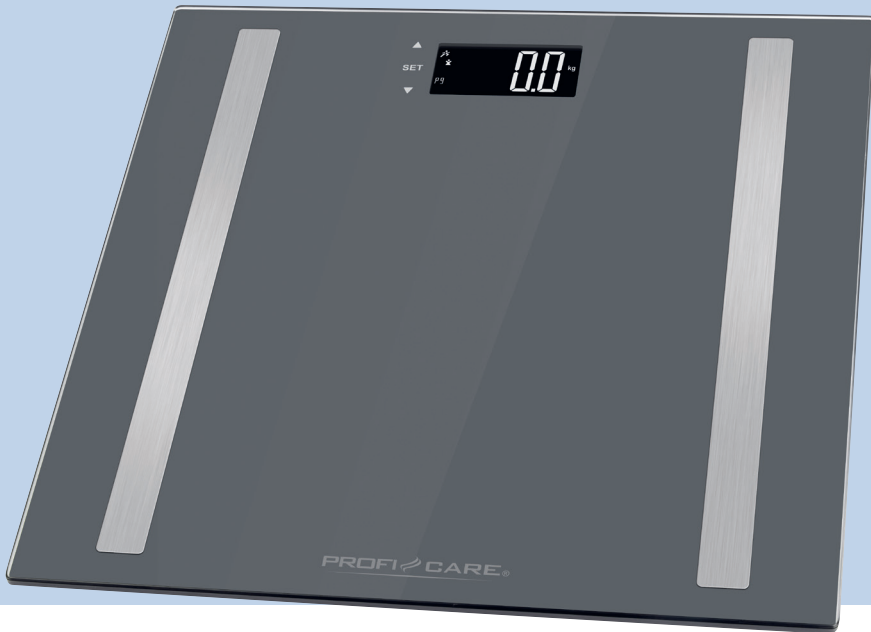## Instruction Manual

Thank you for choosing our product. We hope you will enjoy using the appliance.

### Symbols in this Instruction Manual

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the appliance:

#### WARNING:

This warns you of dangers to your health and indicates possible injury risks.

#### CAUTION:

This refers to possible hazards to the appliance or other objects.

## Content

Overview of the Components.....	3
Display .....	3
Safety Information .....	10
Handling Batteries.....	10
Warnings for the Use of the Appliance .....	11
Unpacking the Appliance.....	11
Overview of the Components /Scope of Delivery .....	11
Inserting the Batteries.....	11
Notes for Use.....	11
Indications in the Display .....	11
Using the Appliance .....	12
Switching On .....	12
Select Unit of Measurement .....	12
Pure Weighing Function.....	12
Saving Personal Data .....	12
Measuring and Analysing.....	12
Switching Off .....	13
Evaluation of the results .....	13
Body Fat (Fat).....	13
Body Water (TBW).....	13
BMI (Body Mass Index).....	13
Muscle Mass (♣).....	14
Bone Weight (♣).....	14
Daily Kilocalorie Requirement (KCAL).....	14
Visceral Fat (♣).....	14
Cleaning.....	14
Storage.....	14
Technical Data.....	15
Disposal .....	15
Battery Disposal.....	15
Meaning of the "Dustbin" Symbol .....	15

## Safety Information

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible, the box with the internal packing. If you pass the appliance on to a third party, always include all the relevant documents.

- The appliance is intended for private use only and not for commercial use.
- Only use the appliance for its intended purpose.
- Keep it away from heat, direct sunlight, moisture (never immerse it in liquids) and sharp edges.
- Do not drop the appliance. This could damage the sensitive electronics.
- Do not repair the appliance by yourself. Always contact an authorized technician.
- Children should be supervised to ensure that they do not play with the appliance.
- In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

#### WARNING:

Do not allow small children to play with the plastic foil.  
**There is a danger of suffocation!**

### Handling Batteries

Observe polarity (+/-)!

Keep away from children!

Non-rechargeable batteries are **not** to be recharged!

Do **not** disassemble batteries!

Do **not** short-circuit batteries!

#### Danger of explosion!

Do **not** expose batteries to intense heat, such as sunshine, fire or similar.

- Batteries can leak battery acid. If the appliance will not be used for a longer period of time, remove the batteries.
- If a battery leaks, do not rub the liquid into the eyes or mucous membranes. In case of contact, wash hands, rinse eyes with clean water. Seek medical attention if any discomfort persists.
- Different battery types and new and used batteries must not be used together.
- Exhausted batteries are to be removed from the appliance and safely disposed of. Batteries are not to be disposed of together with domestic wastes. Return used batteries to the responsible collection points or to the dealer.

## Warnings for the Use of the Appliance

### WARNING:

- If the weighing surface is damp or has drops of water on it, wipe it dry first, as there is a **risk of slipping**.
  - Do not use the scale if you are using a pacemaker or other medical implant with electronic components. The function of the implant may be impaired by the measurement.
- When the scale is being carried or is not in use, please only store it in a horizontal position.
  - The body fat analyser is not suitable for the following groups of people:
    - Children under 10 years of age and adults over 100 years of age
    - Persons with symptoms of fever, oedema or osteoporosis
    - Persons undergoing dialysis treatment
    - People with a pacemaker
    - People taking cardiovascular medication
    - Pregnant women
    - Athletes who do more than 5 hours of intensive exercise per day
    - People with a heart rate of less than 60 beats per minute
    - Persons with significant anatomical deviations in the legs with regard to overall body size (leg length significantly shortened or lengthened).

## Unpacking the Appliance

1. Remove the appliance from its packaging.
2. Remove all packaging material, such as foils, filling material, cable ties and cardboard packaging.
3. Check the appliance for possible transport damage to avoid hazards. Do not operate a damaged appliance. Return it to the dealer immediately.

## Overview of the Components / Scope of Delivery

- 1 UNIT button
- 2 Display
- 3 Contact surfaces
- 4 ▼ Taste
- 5 SET button
- 6 ▲ button
- 7 Battery compartment

## Inserting the Batteries

(Batteries not included)

1. Open the battery compartment on the underside.
2. Insert 2 batteries of type AAA/R03. Please ensure that the polarity is correct (see the floor of the battery compartment)!
  - Do not use rechargeable batteries / accumulators.
3. Close the battery compartment again.

If **Lo** appears in the display, then the batteries must be replaced at once.

## Notes for Use

- The bathroom scale is designed for a maximum weight of 180 kg (396.8 lb).
- To obtain a comparable value, we recommend that you always weigh yourself at the same time of day without clothes and shoes.
- The bathroom scale rounds up or down in steps of 100 g.

## Indications in the Display

All segments in the display appear briefly	Self-test
<b>PO-PG</b>	User memory space
	Male
	Female
<b>cm (ft:in)</b>	Height
<b>age</b>	Age
<b>kg (lb/st:lb)</b>	Body weight
<b>fat</b>	Body fat in %
<b>TBW</b>	Water content in %
<b>BMI</b>	Body mass index
	Muscle mass in %
	Bone weight in kg
<b>KCAL</b>	Basal metabolic rate of the daily energy requirement in kilocalories
	Visceral fat (organ fat)
	Estimated body age
	Sport mode

%	Percentage: body fat, water percentage and muscle mass
Underfat	Underfat
Healthy	Healthy
Overfat	Overfat
Obese	Obese
- - - -	No measuring
UnSt	No stable base
Err2	Exceeding the fat analysis measurement limit or incorrect measurement
Err	Exceeding the weight limit
Lo	Low battery

## Using the Appliance

### Switching On

- Place the bathroom scale on a hard, flat surface (avoid carpeted floors). Uneven floors lead to inaccuracies.
- Weighing function only**  
Switch on the scale by applying light pressure with your foot on the weighing surface.
- Measuring and analysing**  
Press the **SET** button to switch on the scale.
  - All segments are briefly shown on the display.

### Select Unit of Measurement

- Press the **UNIT** button on the back of the scale repeatedly until the desired weight unit is displayed. The unit of measurement for body size changes automatically with.
- The appliance starts in the last unit of measurement selected.

### Pure Weighing Function

- Stand on the scales with both feet. The appliance switches on automatically.
- Distribute your weight evenly and wait while the scale determines your weight.
- The display first flashes and then your exact weight is displayed.
  - When you leave the scale, it switches off automatically after approx. 10 seconds, which extends the battery life.
  - If no weight is determined while "0.0" is displayed, the scale switches off automatically after approx. 10 seconds.

### Saving Personal Data

The scale has preset values that you can change as described below.

Perform the following steps quickly and within a few seconds. Otherwise, the device will automatically exit the input mode.

Press and hold one of the buttons (**▲/▼**) for approx. 2 seconds during the setting to activate fast mode.

- Press the **SET** button. The last selected memory location (e.g. **P0**) flashes on the display. Confirm the setting with **SET** or select a different memory location with the **▲/▼** button. Then press **SET** to move to the next setting.
- The symbol for "male" or "female" flashes. Confirm the setting with **SET** or use the **▲/▼** button to select the symbol for your gender. Also select between normal or sport mode. Then press **SET** to move to the next setting.
- The preset size flashes. Use the **▲** or **▼** button to select your size. Confirm your entry with **SET**.
- The default setting **age** flashes. Select your age using the **▲** or **▼** button. Confirm your entry with **SET**.
- You will then see the memory location, the gender symbol and the weight "0.0 kg" (0.0 lb/0:00 st:lb) on the display.

You now have the following options:

- Press the **SET** button repeatedly to display the settings again or to change the settings.
- You can set additional memory locations. Proceed as described above.
- You can start measuring immediately.

### Measuring and Analysing

The analyser only works correctly if you stand barefoot on the contact surfaces and, if necessary, with slightly moistened soles. Completely dry soles can lead to unsatisfactory results as they have too low conductivity. If necessary, wipe the contact surfaces once to improve the conductive effect.

Save your personal data to a memory location beforehand.

- Press the **SET**. Select the memory location where your data is stored with the **▲/▼** button.
- After selecting the memory location, the stored data is displayed one after the other.
- Stand with both feet on the contact surfaces (3) of the Scale.
- Your weight flashes on the display first. Once the display has stopped flashing, the analysis begins. During the analysis you will see a treadmill "0000" on the display. Remain on the scales until the last display (**♂**) is displayed.

The following results are displayed in sequence:

- The weight in the previously selected unit (**kg/lb/st:lb**)
- The body fat percentage (**fat**) + weight status display
- The body water percentage (**TBW**)
- The body mass index (**BMI**)

- 5 The muscle mass percentage (♣)
- 6 The bone weight in kg (♣)
- 7 The basal metabolic rate of the daily energy requirement in **KCAL**
- 8 Visceral fat (♣)
- 9 The estimated body age (♣)

The values determined are then repeated and the scales switch off automatically.

### Switching Off

The automatic switch-off function switches the scale off automatically after approx. 10 seconds.

## Evaluation of the results

Our body consists of body fat, body water, muscle mass and bone weight, among other things. Due to individual water retention in these main components, the values determined may vary. This is normal and strongly dependent on diet and lifestyle.

There are heavy people who have a balanced ratio between fat and muscle mass and there are apparently slim people with significantly increased body fat values. Fat is a factor that contributes significantly to the strain on the organism, especially the cardiovascular system. Controlling both values (fat and weight) is therefore a basic prerequisite for health and fitness.

The calculated measurement results of body fat and body water percentage are only reference values. They are not recommended for use for medical purposes. Please always consult your doctor if you have any questions about medical treatment or a diet to achieve a normal weight.

### Body Fat (Fat)

Body fat is an important component of the organism. From protecting the joints and storing vitamins to regulating body temperature, fat plays an important role in the body. The aim is therefore not to radically reduce your body fat percentage, but to achieve and maintain a balanced ratio between muscle mass and fat mass.

You can find reference values for your fitness profile in the guide value table. It shows the dependence of body fat values on age and pointer.

*Guide values for body fat percentage:*

Age	Women			
	Underfat	Healthy	Overfat	Obese
18–39	≤ 20.9	21–34.9	35–39.9	≥ 40
40–59	≤ 21.9	22–35.9	36–40.9	≥ 41
> 60	≤ 22.9	23–36.9	37–41.9	≥ 42

Age	Men			
	Underfat	Healthy	Overfat	Obese
18–39	≤ 10.9	11–21.9	22–26.9	≥ 27
40–59	≤ 11.9	12–22.9	23–27.9	≥ 28
> 60	≤ 13.9	14–24.9	25–29.9	≥ 30

### Body Water (TBW)

The human organism consists of around 55–60% water, in varying amounts depending on age and gender.

Water has a number of tasks to fulfil in the body:

- It is a building block of our cells. This means that all cells in the body, whether skin, glandular, muscle, brain cells or others, can only function properly if they contain enough water.
- It is a solvent for important substances in our organism.
- It is a means of transport for nutrients, endogenous substances and metabolic products.

*Guide values for body water content in per cent:*

Age	Women		
	Underfat	Healthy	Overfat
1–30	< 49.5	49.6–52.9	> 53
> 31	< 48.1	48.2–51.5	> 51.6
Age	Men		
	Underfat	Healthy	Overfat
1–30	< 53.6	53.7–57	> 57.1
> 31	< 52.3	52.4–55.6	> 55.7

### BMI (Body Mass Index)

The Body Mass Index (BMI) is the ratio of body weight in kilograms and height in metres squared.

BMI = weight in kg : (height in metres x height in metres)

*Example:*

A woman weighs 55 kg and is 1.68 metres tall. Her BMI is calculated as follows:

55 kg : (1.68 m x 1.68 m) ≈ 19.5 kg/m<sup>2</sup> ⇨ normal weight

The "height squared" value is not related to body surface area. The BMI is only a rough guide, as it does not take into account age and gender or the individual composition of a person's body mass of fat and muscle tissue.

*Guide values for BMI:*

BMI	Underfat	Healthy	Overfat	Obese
		< 18.5	18.6–24	24.1–28

Note that for very muscular bodies (body builders), the interpretation of the BMI indicates overweight. The reason for this is that muscle mass, which is well above average, is not carried into account in the BMI formula.

### Muscle Mass (👤)

The muscle mass percentage is the percentage of muscle mass in relation to total body weight. The ideal muscle mass percentage depends on gender, age and training condition.

*Guide values for muscle mass percentage:*

Women		
Underfat	Healthy	Overfat
< 25	25.1–30	> 30.1
Men		
Underfat	Healthy	Overfat
< 31	31.1–39	> 39.1

### Bone Weight (👤)

Bone weight refers to the total weight of all bones in the body. This weight varies depending on gender, age and body size.

*Guideline values for bone weight in kilograms:*

Women			
Age	Underfat	Healthy	Overfat
1–39	< 1.7	1.8–2.2	> 2.3
40–60	< 2.1	2.2–2.7	> 2.8
> 61	< 2.4	2.5–3.1	> 3.2
Men			
Age	Underfat	Healthy	Overfat
1–54	< 2.4	2.5–3.1	> 2.8
55–75	< 2.8	2.9–3.6	> 3.7
> 76	< 3.1	3.2–4.0	> 4.1

### Daily Kilocalorie Requirement (KCAL)

The daily kilocalorie requirement depends on age, gender, weight and height. It is made up of the basal metabolic rate (energy at rest) and the power metabolic rate (energy for activities).

The basal metabolic rate is the amount of energy that the body needs at complete rest for vital functions such as breathing, heartbeat and temperature regulation.

The power metabolic rate is the additional energy that the body uses for physical and mental activities.

*Guideline values for daily kilocalorie requirements:*

	Weight	Age			
		< 17	18–29	30–50	> 50
Women	< 50 kg	1368	1311	1213	1120
	50–65 kg	1511	1455	1357	1263
	65–80 kg	1654	1598	1457	1407
	80–95 kg	1798	1742	1643	1550
	> 95 kg	2180	2124	2026	1836
Men	Weight	Age			
		< 17	18–29	30–50	> 50
	< 50 kg	1489	1408	1266	1131
	50–65 kg	1695	1614	1473	1337
	65–80 kg	1902	1821	1679	1544
80–95 kg	2108	2027	1885	1750	
> 95 kg	2864	2783	2435	2300	

### Visceral Fat (👤)

Visceral fat refers to the fat stored in the free abdominal cavity that surrounds the internal organs. It protects the organs and serves as an energy store.

*Guideline values for visceral fat percentage:*

Normal	A little too much	Too much
< 9	10–14	> 15

## Cleaning

### ⚠️ CAUTION:

- Do **not** immerse the appliance in water! This could destroy the electronics.
- Do not use a wire brush or other abrasive objects.
- Do not use any sharp or abrasive cleaning agents.
- Clean the appliance with a slightly damp, not wet, cloth. Make sure that no moisture or liquid enters the appliance!
- You can use a commercially available washing-up liquid or a simple soap as an additive.
- After cleaning, dry the appliance well with a dry, soft cloth.

## Storage

Store the appliance in a clean, dry place.

## Technical Data

Model:..... PC-PW 3149 FA  
 Net weight: ..... approx. 1.54 kg  
 Battery: ..... 2 × 1.5 V Type AAA/R03  
 Memory locations:..... 10  
 Height adjustment: .... from 80 to 240 cm (27 ½"–7'10.5" ft.in)  
 Age: ..... from 10 to 100 years  
 Minimum weight: ..... from 3 kg (6.6 lb/0: 9.8 st.lb)  
 Maximum weight: ..... 180 kg (396.8 lb/28: 3 ½ st.lb)  
 Resolution for  
 weight measurement: ..... 100 g (0.2 lb/¼ st.lb) steps  
 Body fat graduation: ..... 0,1 %  
 Body water graduation: ..... 0,1 %  
 Body fat measurement range: ..... 1 %–60 %  
 Body water measurement range: ..... 20 %–85 %  
 Muscle mass measurement range: ..... 10 %–95 %  
 Bone weight measurement range: ..... 0.5 kg–10 kg

The right to make technical and design modifications in the course of continuous product development remains reserved.

This appliance has been tested according to all applicable, current CE directives and built according to the latest safety regulations.

dealer and contractual partner is also obliged to take back the old appliance free of charge.

## Disposal

### Battery Disposal

As a consumer you are obliged by law to return old batteries / accumulators.

Your appliance contains batteries. If these batteries are spent, follow the instructions below:

#### CAUTION:

Batteries / rechargeable batteries must not be disposed of with household waste!

- The batteries must be removed from the appliance.
- Take spent batteries to the appropriate collection point or to a dealer.
- Your town or local authority can provide information about public collection points.

### Meaning of the "Dustbin" Symbol

Protect our environment: do not dispose of electrical equipment in the domestic waste.

Please return any electrical equipment that you will no longer use to the collection points provided for their disposal.

This helps avoid the potential effects of incorrect disposal on the environment and human health.

This will contribute to the recycling and other forms of re-utilisation of electrical and electronic equipment.

Information concerning where the equipment can be disposed of can be obtained from your local authority. Your

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